

Isolation and identification of restriction endonuclease SseAI

M.Rina, A.Karagouni¹, M.Pagomenou, I.Tsigos and V.Bouriotis*

Institute of Molecular Biology and Biotechnology, Enzyme Technology Division, PO Box 1515, Heraklion 711 10, Crete and ¹Institute of General Botany, University of Athens, Athens 15784, Greece

Submitted July 5, 1991

SseAI, an isoschizomer of NarI (1), has been purified from *Streptomyces* species. SseAI recognises the palindromic sequence 5'...GGCGCC...3' and cleaves between second G and C. The enzyme was purified using the following chromatographic steps. 1. Phosphocellulose, 2. Heparin-Sepharose. The enzyme was free of contaminating nuclease activity. After 100 fold overdigestion on Adeno-2 DNA greater than 95% of the DNA fragments can be ligated and recut by SseAI. Optimal conditions for enzyme activity are 6 mM Tris-HCl (pH 7.4), 6 mM MgCl₂, 6 mM β-mercaptoethanol at 37°C. The fragments produced by SseAI digestion of lambda-HindIII DNA, Adeno-2 DNA, SV40, Φx174, pBR322 and lambda DNA match those predicted by cleavage at the sequence GGCGCC (figure 1, lanes 4–9).

In order to determine the cleavage site within the recognition sequence a pUC19 vector which contained a recognition site for the enzyme was digested by the enzyme SseAI then annealed with sequencing primers and extended with Klenow enzyme in the presence of α³²P-dATP. Dideoxy sequencing reactions were performed at this region with the same primers and run in parallel with the extended products (2). Results in figure 2 show that the extended products of both primers (lanes F and R) comigrate with the band corresponding to the 3'G in the 5'...GGCGCC...3' sequence. Therefore SseAI recognises and cleaves the following sequence.



REFERENCES

1. Comp,D., Wilson,G., Schildkraut,I. and Greenough,L. Cited in Roberts,R.J. (1985) *Nucl. Acids Res.* **13**, 165–200.
2. Sanger,F., Nicklen,S. and Coulson,A.R. (1977) *Proc. Natl. Acad. Sci. USA* **74**, 5463–5467.

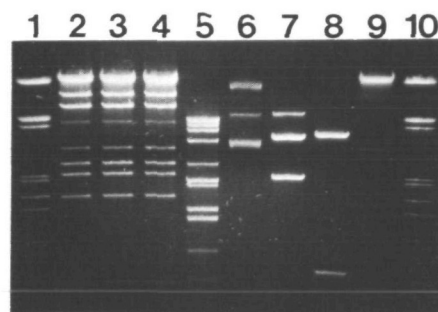


Figure 1. SseAI digests: lane 2: lambda-HindIII DNA digested by NarI, 3: lambda-HindIII DNA digested by NarI and SseAI, 4: lambda-HindIII DNA, 5: Adeno-2, 6: SV40, 7: Φx174, 8: pBR322, 9: lambda DNA, 1, 10: lambda HindIII-EcoRI size standard.

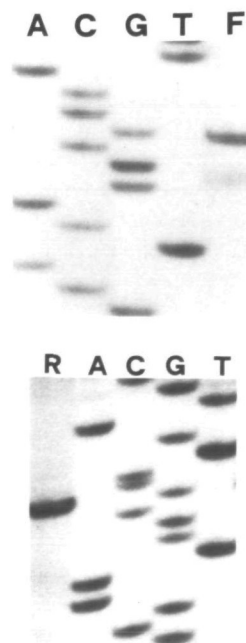


Figure 2.